
ZEEN CARE

BOOST PRO MASSAGE GUN

ZC-JM05

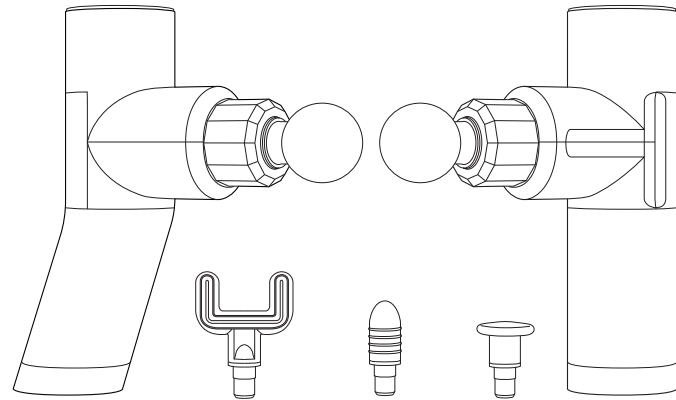Upute za korištenje

Boost Pro massage gun

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Cleaning and Maintenance

1. Always unplug the charger from the power outlet before cleaning.
2. Wipe the device with a neutral detergent and a soft cloth. The use of corrosive cleaning agents is strictly prohibited.
3. Prevent water or other corrosive liquids from entering the device, as this may cause malfunction or damage.

Safety tips

Do not open or use the product if it is damaged. Feniks Info d.o.o. accepts no responsibility for product damage or injuries caused by careless, improper, or unintended use, or by use not in accordance with the intended purpose of the product. Use the device only in closed and dry indoor areas. There is a risk of short circuit if it comes into contact with moisture. It is intended exclusively for private use in accordance with these instructions and must not be used for commercial purposes.

DO NOT open or use the device if it is damaged!

DO NOT use it in damp or wet locations!

Only charge the device within the specified technical parameters!

DO NOT use the device at high temperatures and/or if it is covered and/or in a humid environment!

DO NOT expose the device to direct sunlight!

DO NOT overload the device!

Always unplug the charger from the outlet before cleaning, maintenance, or using the device!

If any malfunction occurs (e.g. burning smell), immediately turn off the device and stop using it!

DO NOT touch the charger with wet hands!

DO NOT use the device in the bathroom or other humid areas!

DO NOT spray volatile oils, coatings, or pesticides on the device!

DO NOT use damaged chargers or cables!

Keep children away while the device is in operation!

Device Disposal

This device must not be disposed of with household waste. Please take it to a designated collection point for electronic waste in your community or district.

Packaging Disposal

The packaging must not be disposed of with household waste. Please follow the appropriate local procedures for disposal.

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Guide for use

Boost Pro massage gun


Please read this manual carefully BEFORE using the device to ensure safe and reliable operation. Keep the instructions in a safe place and pass them on to any future owners.

Technical Specifications

Product Name	Boost Pro Massage Gun
Model	ZC-JM05
Input	USB C / 5 V / 24 W / 1A
USB Cable	0,5 m
Speeds	4
Manual Operation	Yes
Program Duration	15 minutes
Rotation Speed	4000 ± 200 rpm
Number of Attachments	4
Charging Time	2,5 hours
Amplitude	0,6 cm
Battery	18650, 7,4 V, 1200 mAh
Battery Life	260~350 minutes
Warranty	12 months

Instructions for Use

1. Press and hold the button on the device to turn it on. Once turned on, the device will automatically be set to the first speed, and one red indicator LED will light up.
2. Press the button again to switch to the second speed, and two blue indicator LEDs will light up.
3. With each additional press of the button, the device switches to a higher speed, and the number of illuminated LEDs increases accordingly.
4. The device has a total of four speeds. When it reaches the fourth speed, pressing the button again will return it to the first speed.
5. Press and hold the button on the device to turn it off.
5. Pritisnite i držite tipku na uređaju kako biste isključili uređaj.

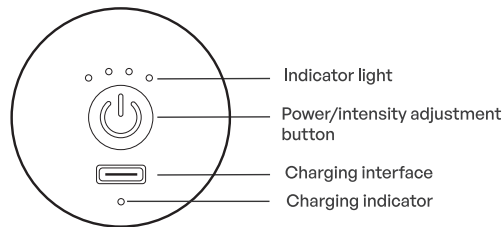
Charging Instructions

1. Do not use the device while it is connected to the charger.
2. Connect the USB-C connector of the included USB cable to the USB-C port on the massager.
3. Connect the USB-A plug of the included USB cable to a 5V, minimum 1A USB charger (charger not included).
4. Charging will start automatically and the red LED indicator will light up while charging.
5. Once the device is fully charged, the LED will turn green. When charging is complete, you can disconnect the charger. Disconnect in the reverse order of connection.

Method of replacement of massage head

1. Pull out the original massage head.

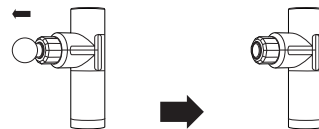2. Put on the massage head you need.