

Menjava masažnih glav

1. Odstranite obstoječo masažno glavo.

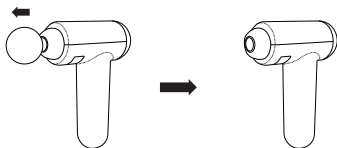2. Namestite novo glavo, ki jo želite uporabiti.

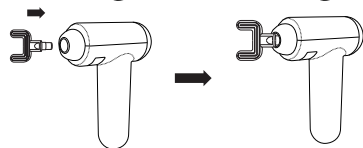ZEEN CARE

PULSE MASSAGE GUN

ZC-JM03

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

Method of replacement of massage head

1. Pull out the original massage head.

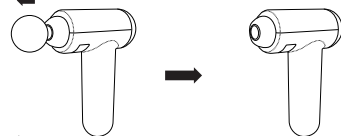2. Put on the massage head you need.

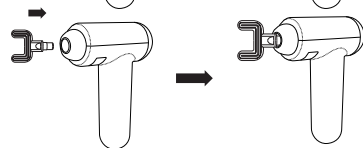⚠ Safety tips

Do not open or use the product if it is damaged. Feniks Info d.o.o. accepts no responsibility for any damage to the product or injuries resulting from careless, improper, or unintended use, or use that is not in accordance with the intended purpose of the product. Use the device only in indoor and dry environments. There is a risk of short circuit if it comes into contact with moisture. This device is intended exclusively for private use in accordance with these instructions and must not be used for commercial purposes.

DO NOT open or use the device if it is damaged!

DO NOT use it in damp or wet areas!

Charge the device only within the specified technical parameters!

DO NOT use the device at high temperatures and/or if it is covered and/or in humid environments!

DO NOT expose the device to direct sunlight!

DO NOT overload the device!

Always unplug the charger from the socket before cleaning, maintenance, or using the device!

In case of malfunction (e.g., burning smell), immediately switch off the device and stop using it!

Do not touch the charger with wet hands!

Do not use the device in the bathroom or other humid rooms!

Do not spray volatile oils, coatings, or pesticides on the device!

Do not use damaged chargers or cables!

Keep children away from the device while it is operating!

Device Disposal - This device must not be disposed of with household waste. Please take it to a designated collection point for electronic waste in your community or district.

Packaging Disposal - The packaging must not be disposed of with household waste. Please follow the appropriate local procedures for disposal.

Massage heads and Feature

Round head

Used for hand, back, buttock, leg massage.

Flat head

Used for relaxing and shaping each muscle part.

Taper head

Used for deep tissue, such as joints, palm and plantar.

U-shaped head

Used for neck, spinal and achilles massage.

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
Instructions for Use

1. Press the button on the device to turn it on. After switching on, the device will automatically start at the first speed and one LED indicator will light up red.
2. Press the button again to switch to the second speed, and two LED indicators will light up red.
3. Each subsequent press of the button increases the speed, and the number of lit LED indicators increases accordingly.
4. The device has a total of six speeds. When it reaches the sixth speed, press the button again to turn it off.
5. You can turn off the device at any time by pressing and holding the button.

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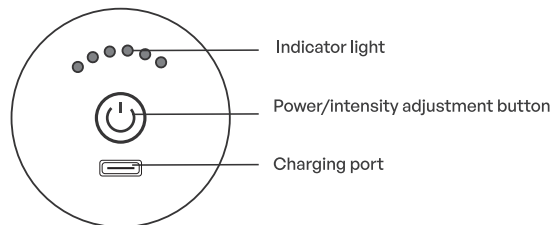Guide for use

Pulse massage gun

Please read this manual carefully BEFORE using the device to ensure safe and reliable operation. Keep the instructions in a safe place and pass them on to any future owners.

Technical Specifications

Product Name	Pulse Massage Gun
Model	ZC-JM03
Input	USB C / 5 V / 24 W / 1A
USB Cable	0,5 m
Speeds	6
Manual Operation	Yes
Program Duration	15 minutes
Rotation Speed	3000 ± 200 rpm
Number of Attachments	4
Charging Time	2,5 hours
Amplitude	0,6 cm
Battery	18650, 7,4 V, 1200 mAh
Battery Life	150–180 minutes
Warranty	12 months

Charging Instructions

1. Do not use the device while it is connected to the charger.
2. Connect the USB C connector of the included USB cable to the USB C input on the massager.
3. Connect the USB A plug of the included USB cable to a 5V, minimum 1A USB charger (charger not included).
4. Charging will start automatically, and the red LED light will turn on while charging.
5. When the device is fully charged, the red LED light will turn off. Once charging is complete, you may disconnect the charger. Disconnect in reverse order.

Cleaning and Maintenance

1. Always unplug the charger from the socket before cleaning.
2. Wipe the device with a neutral detergent and a soft cloth. Do not use corrosive cleaning agents.
3. Prevent water or other corrosive liquids from entering the device, as this may cause malfunction or damage.